**Bestand<sup>®</sup>**

**PERCUSSIVE THERAPY DEVICE**

## Feature Guide

High-power motor with a maximum thrust reaching 13kgs; the first choice for a sportsman.


Ergonomic handle design, comfortable to hold and use.

Visualized gear display, for a more convenient use.

Innovative transmission mechanism, lower noise than traditional products.


Optional massage heads, which allow every part of your body to experience relaxation.

Large-capacity batteries, longer service lifespan.
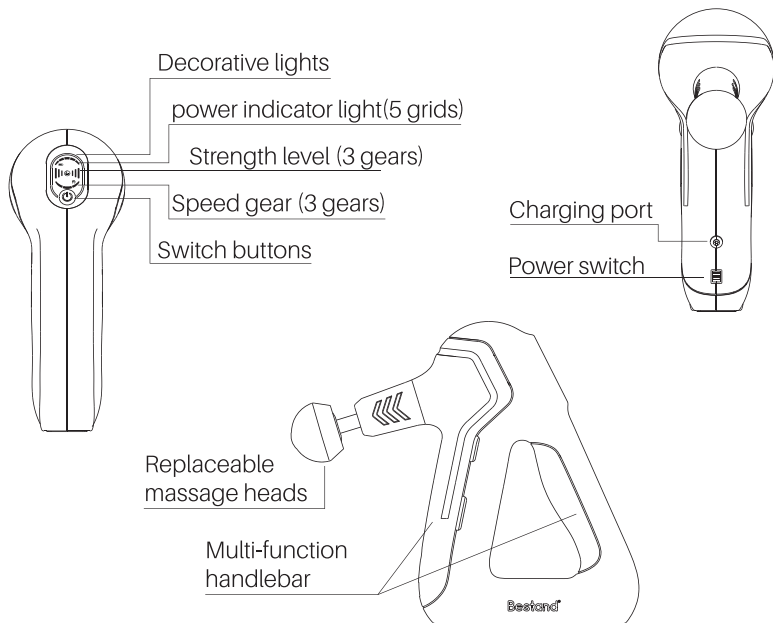
## How to use

1. Turn the power switch to the "ON" position, to turn the power on; at the same time, the decorative light (which flashes once) and the power indicator on the display panel will turn on.

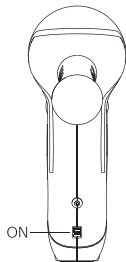2. Pressing the gear switch button will start the motor, and the massage head will start to run in 1st gear; press again to run in 2nd gear; press 3 times to run in 3rd gear; press 4 times to turn off the motor, which causes the massage head to stop running, and places the product in a standby state.


3. After using the product, please turn the power switch to the "OFF" position, which turns the power off, as well as turning all indicators off.

When charging, first insert the charging plug of the charger into the charging port, and then connect the charger to the power socket; when power indicator flashes, the product is charging normally.


After fully charging, the power indicators will all be on.

Note: When charging, the power remains off and the product cannot be used normally.


## Massage head

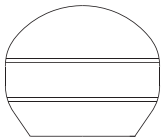**Large sphere**  
suitable for large muscle  
groups

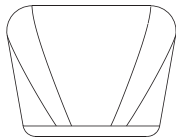**Wedge-shaped head**  
scapular muscle group or  
iliotibial band

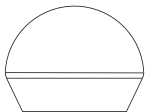**Medium sphere**  
comprehensive massage

**Finger-shaped head**  
(included)  
trigger point massage

Connecting the massage head: make sure the massage head is aligned with the head of the exercise arm, and push it to connect. To remove it, just pull out the massage head.

## Specifications

Product model	MC06	5 grid power indicator	10%-25%-50%-75%-100%
Product weight	1150g	3 gear vibration speed	50%-75%-100%
Charging voltage and current	12.8V, 1.5A	Package contents	massager, massage head, charger, manual
Output power	40W		
Motor power	72W		
Battery capacity	2500mAh		
Charging time	4-4.5 hours		
Battery life	1-6 hours		
3 gear massage intensity	<p>Idling the motor or heavily pressing the telescopic head may cause the intensity indicator to flash. When applying pressure with the telescopic head, the strength indicator remains on all the time.</p>		

# 12 months warranty

Within 12 months of purchasing this product, if you find any problems during use, please feel free to contact us.

Official website: [www.bestand.com](http://www.bestand.com)

Email: [globalbiz@bestand.com](mailto:globalbiz@bestand.com)

## Tips:

If the product cannot be used, it might be due to the battery running out of charge. Please ensure the product is charged before use; when not in use, turn off the power switch.

Please clean and disinfect the product regularly to avoid bacterial growth. Cleaning liquid is recommended to be used for disinfection. Use a cotton cloth, dust-free paper towels, etc. to wipe and allow to dry naturally, avoiding exposure to the sun.

If the product is damaged or has sharp edges before it has been used, please stop using this product.

If you feel unwell during use, please stop using the product immediately to avoid hurting yourself.

This product is a health-care massager for adults, not a toy for children; please store it properly and keep it out of reach of children and pets; if this product or its accessories are swallowed by children or pets, please go to a nearby hospital immediately.

## Warning:

1. This product is not water-resistant, please do not immerse it in water or wash it with water; just wipe it with a wrung-out wet towel. Otherwise, it may cause a short circuit, thus damaging the product.
2. Please take care when using this product, as it is equipped with high-quality polymer batteries. Avoid bashing the product, do not bend it or puncture any part of it by yourself, do not expose this product to fire or high temperatures, and if the product is damaged, bulging, or in any other adverse condition, please stop using it.
3. During operation, please do not touch or insert your fingers or other objects into the gap between the massage head and the body of the product, to avoid pinching.
4. Do not leave the charger plugged in for a long time, since it may cause leakage or other serious accidents due to ageing and other reasons.
5. Do not use the product on the same part of your body for a long time.
6. The power of this product is relatively high. It is recommended for use for only 10-20 minutes at a time to avoid physical discomfort.
7. Do not use this product on fragile bones and weak muscles.
8. Infants, minors, pregnant women, those who wear a pacemaker, have broken bones, or have metal contents in their bodies, please do not use this product.
9. If you feel unwell, stop using it immediately.
10. This product is subject to change without notice.

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